



## FiltaNeb™

Cirrus™2 nebuliser, adult, Intersurgical EcoLite™ non-vented mask kit with filter and tube



Oxygen and Aerosol Therapy ▪ Nebuliser Therapy

To meet the changing requirements of the global healthcare market this product option has been developed to reduce particulates entering the environment during nebuliser therapy. The FiltaNeb™ features a non-vented mask based upon the Intersurgical EcoLite™ mask, with a filter on the exhalation gas pathway. The mask seal is coloured green to distinguish this non-vented mask from standard aerosol masks. To ensure that the best possible fit can be achieved, the mask kit includes an additional head strap.

The system provides >95% reduction of particulates into the environment, helping to minimise potentially infected exhaled gases and emissions escaping during nebulisation therapy, reducing the risk to healthcare professionals.

The Cirrus™2 nebuliser incorporates a number of features, including a graduated drug cup that provides guidance on how much drug is in the nebuliser. The cup has a quick-release fitting which enables the nebuliser to be attached, as well as detached, with a quarter turn. This speeds up the filling process and provides a more reliable seal between the cup and the nebuliser body.

The Cirrus2 nebuliser is designed to deliver drugs for tracheobronchial deposition.

**Incurved nose seal**

Conforms to different nose shapes designed to reduce the risk of aerosolised drug entering patient's eyes

**Filter on exhalation port**

Reduces particulates entering the environment during nebuliser therapy

**Graduated nebuliser cup guide**

Faster preparation

**Low residual volume**

Reduced drug waste

**Additional head strap**

provides an improved seal between the mask and the patient's face

**Soft face seals**

Increased patient comfort

**Green mask seal**


Differentiates the non-vented mask from standard aerosol masks

**A choice of "under chin" or "on chin" positions**

Provides a better fit on a wider range of patient face shapes

**Quick-release fitting**

Attaches and detaches easily

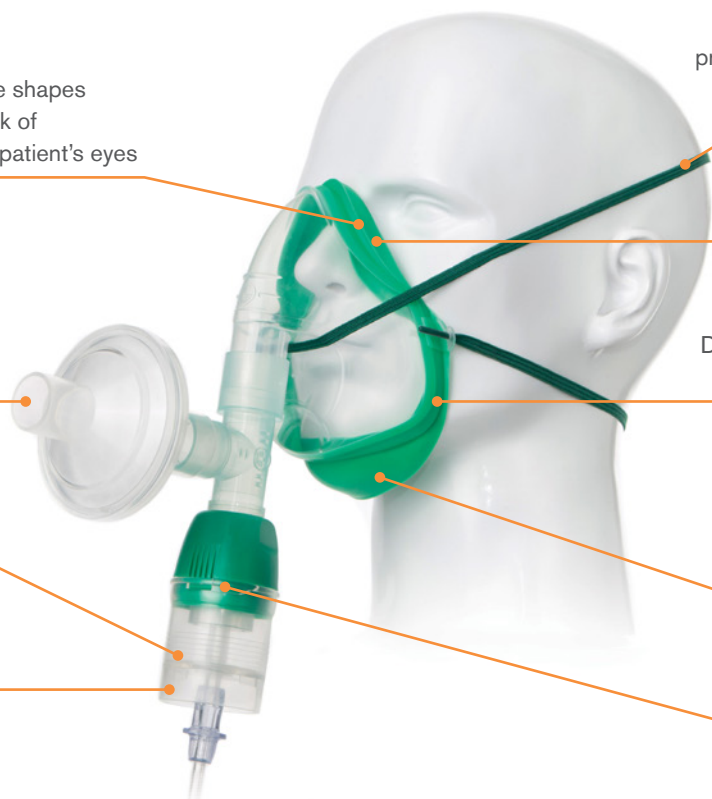

Code	Description	Tube length	Box Qty.
1453044	FiltaNeb™ Cirrus™2 nebuliser, adult, Intersurgical EcoLite™ non-vented mask kit with filter and tube	2.1m	12

Lower environmental impact product

IS10.39 • Issue 4 11.23



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

**Please think before you print**  
Save energy and paper.  
If you must print this information sheet please print it double sided.