




# O2-MAX™

Disposable CPAP device

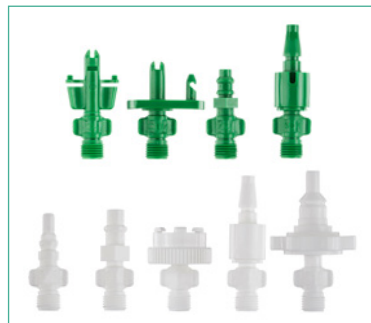

Critical Care ▪ Pulmodyne®



Quality, innovation and trust

Interact with us




[www.intersurgical.com](http://www.intersurgical.com)

## O2-MAX™

The O2-MAX System can deliver up to 140 L/min of flow to the patient, with 30% FiO<sub>2</sub>. Customise your FiO<sub>2</sub> to increase from 30, 60 or 90% with the Trio System. Adjustable PEEP options includes a 3-SET™ (5/7.5/10 cm H<sub>2</sub>O) or 5-SET™ (5/7.5/10/12.5/15 cm H<sub>2</sub>O) to suit a variety of patient needs.

Additionally, this completely disposable system offers integrated nebulisation while maintaining a consistent flow to the patient. A variety of Quik-Connect™ allow direct fitting of the O2-MAX to wall oxygen for immediate therapy activation.

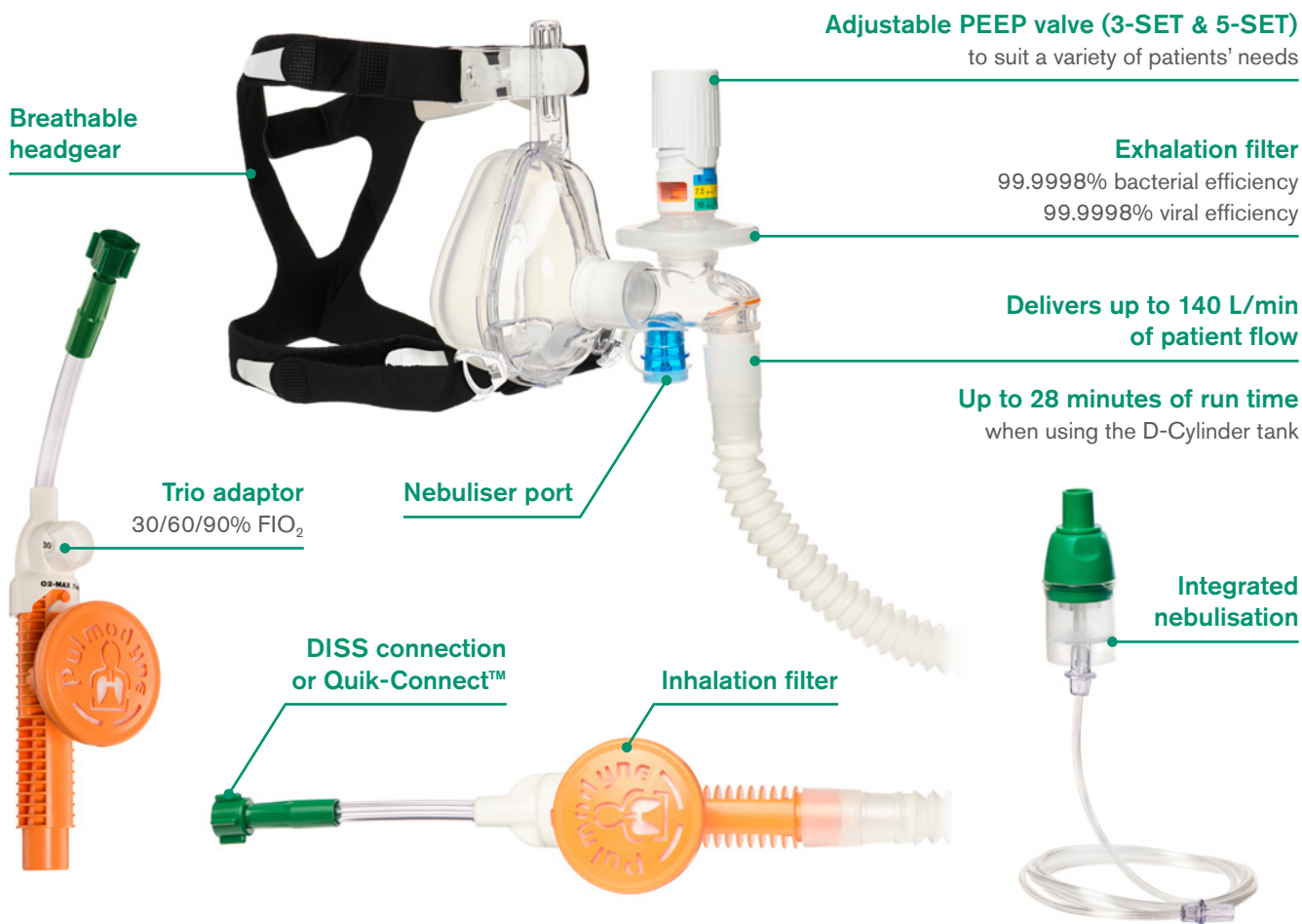

### O2-MAX accessories

Code	Description	Box Qty.
313-8000	O2-MAX Trio control unit	10
313-8300	O2-MAX exhalation filter	10
313-7716	O2-CPAP™ fixed value PEEP valve, 5.0 cm H <sub>2</sub> O, blue	10
313-7717	O2-CPAP fixed value PEEP valve, 7.5 cm H <sub>2</sub> O, yellow	10
313-7718	O2-CPAP fixed value PEEP valve, 10.0 cm H <sub>2</sub> O, green	10
313-7719	O2-CPAP fixed value PEEP valve, 12.5 cm H <sub>2</sub> O, pink	10
313-7720	O2-CPAP fixed value PEEP valve, 15.0 cm H <sub>2</sub> O, red	10
313-7721	O2-CPAP fixed value PEEP valve, 20.0 cm H <sub>2</sub> O, gold	10
313-7725	3-SET adjustable PEEP valve (5.0/7.5/10.0 cm H <sub>2</sub> O), white	10
313-7730	5-SET adjustable PEEP valve (5.0/7.5/10.0/12.5/15.0 cm H <sub>2</sub> O), grey	10

Make an enquiry



## O2-MAX™ Fixed System with expandable tube (delivers 30% FiO<sub>2</sub>)

[Make an enquiry](#)

Code	Description	Box Qty.
313-7555XN	O2-MAX Fixed System with expandable tube, 3-SET adjustable PEEP valve, nebuliser and BiTrac™ ED mask, large adult	10
313-7556XN	O2-MAX Fixed System with expandable tube, 3-SET adjustable PEEP valve, nebuliser and BiTrac ED mask, medium adult	10
313-7555X	O2-MAX Fixed System with expandable tube, 3-SET adjustable PEEP valve and BiTrac ED mask, large adult	10
313-7556X	O2-MAX Fixed System with expandable tube, 3-SET adjustable PEEP valve and BiTrac ED mask, medium adult	10
313-7554XN	O2-MAX Fixed System with expandable tube, 5-SET adjustable PEEP valve, nebuliser and BiTrac ED mask, large adult	10
313-7557XN	O2-MAX Fixed System with expandable tube, 5-SET adjustable PEEP valve, nebuliser and BiTrac ED mask, medium adult	10
313-7554X	O2-MAX Fixed System with expandable tube, 5-SET adjustable PEEP valve and BiTrac ED mask, large adult	10
313-7557X	O2-MAX Fixed System with expandable tube, 5-SET adjustable PEEP valve and BiTrac ED mask, medium adult	10

## O2-MAX Trio System (delivers adjustable 30%, 60%, 90% FiO<sub>2</sub>)

[Make an enquiry](#)

Code	Description	Box Qty.
313-8018XN	O2-MAX Trio System with expandable tube, 5-SET adjustable PEEP valve, nebuliser and BiTrac ED mask, large adult	10
313-8008XN	O2-MAX Trio System with expandable tube, 5-SET adjustable PEEP valve, nebuliser and BiTrac ED mask, medium adult	10
313-8018	O2-MAX Trio System with smoothbore tube, 5-SET adjustable PEEP valve and BiTrac ED mask, large adult	10
313-8008	O2-MAX Trio System with smoothbore tube, 5-SET adjustable PEEP valve and BiTrac ED mask, medium adult	10

## Quik-Connect™ disposable oxygen connections

[Make an enquiry](#)


Code	Description	Box Qty.
313-111111	Quik-Connect Ohmeda gas fitting, green	30
313-121111	Quik-Connect Chemetron gas fitting, green	30
313-131111	Quik-Connect Puritan Bennett gas fitting, green	30
313-141111	Quik-Connect Schrader gas fitting, green	30
313-351111	Quik-Connect AGA gas fitting, white	30
313-361111	Quik-Connect DIN gas fitting, white	30
313-371111	Quik-Connect AFNOR gas fitting, white	30
313-381111	Quik-Connect BS gas fitting, white	30
313-391111	Quik-Connect UNI gas fitting, white	30

## O2-MAX™ accessories

The BiTrac ED mask can be supplied separately to integrate the patient interface size range for O2-MAX. This mask features a dual lip seal and the adjustable Omni-Clip™ to guarantee optimal fit. The connection with the breathing system is a 30M/22F port that fits directly to the O2-MAX patient connections as well as to the ADAPT-2-VENT™ elbows accessories.



### BiTrac™ ED mask with CAP headgear

Code	Description	Box Qty.
313-7030C	BiTrac ED mask with CAP headgear, large adult	10
313-7029C	BiTrac ED mask with CAP headgear, medium adult	10
313-7028C	BiTrac ED mask with CAP headgear, small adult	10

[Make an enquiry](#)



### BiTrac ED mask with standard headgear

Code	Description	Box Qty.
313-7031	BiTrac ED mask, extra large adult	10
313-7030	BiTrac ED mask, large adult	10
313-7029	BiTrac ED mask, medium adult	10
313-7028	BiTrac ED mask, small adult	10

[Make an enquiry](#)

### ADAPT-2-VENT™ and Cirrus™2 nebuliser

The ADAPT-2-VENT is an adapter elbow that allows the use of the BiTrac ED mask on standard dual limb or single limb ventilator circuits. This allows the patient to move from CPAP, using the O2-MAX to other forms of ventilation without changing the patient interface. The 18M nebuliser that fits the O2-MAX can be supplied separately to deliver nebulised drugs during CPAP therapy.



313-9185



313-9186

Code	Description	Box Qty.
313-9185	ADAPT-2-VENT 30 mm F non-vented elbow with anti-asphyxiation valve	10
313-9186	ADAPT-2-VENT 30 mm F standard non-vented elbow	10
1402018	Cirrus™2 nebuliser 18M with tube, 1.8 m	45

3-SET, 5-SET, ADAPT-2-VENT, BiTrac, O2-CPAP, O2-MAX and Quik-Connect are manufactured by Pulmodyne Inc and distributed by Intersurgical Ltd.

Ohmeda is a registered trademark of GE Healthcare Finland Oy.

Puritan Bennett is a registered trademark of Covidien LP.

Chemetron is a registered trademark of Allied Healthcare Products, Inc.

**Note:** Not all products are available in all markets, please contact your local office for details.

IS1.64 ▪ Issue 5 02.26



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com

UK ▪ Australia ▪ Benelux ▪ Brazil ▪ Canada ▪ China ▪ Colombia ▪ Czech Republic ▪ Denmark ▪ France ▪ Germany ▪ India ▪ Ireland ▪ Italy  
Japan ▪ Lithuania ▪ Malaysia ▪ Mexico ▪ Norway ▪ Philippines ▪ Portugal ▪ Russia ▪ South Africa ▪ Spain ▪ Sweden ▪ Taiwan ▪ Turkey ▪ USA