



## HOT Top™2

A range of high output nebulisers to reduce treatment times



Oxygen and Aerosol Therapy ▪ Nebuliser Therapy



## HOT Top™2 nebuliser

The HOT Top2 nebuliser is an open nebuliser incorporating twin entrainment ports for high output flows resulting in reduced treatment times which may improve patient compliance.

In addition, the HOT Top2 nebuliser includes a graduated nebuliser cup which provides guidance on how much medication is left in the nebuliser and a quick-release fitting which enables the nebuliser to be attached and detached with a quarter-turn. This speeds up the refilling process and provides a reliable seal between the cup and the nebuliser body.

The HOT Top2 nebuliser is designed to deliver medication for tracheobronchial deposition. Nebuliser performance specifications are detailed in documentation provided with the products, and are available upon request.

At the lower driving gas flow of 6 L/min the performance is similar to a conventional closed nebuliser, making it an option for use with devices which have lower flow outputs.

The HOT Top2 mask kits incorporate the Intersurgical EcoLite™ aerosol mask, this non-PVC mask reduces the environmental impact of this product.

[▶ Watch the video](#)

### Standard 22F taper

Fits masks, T-pieces and mouthpieces

### Quick-release fitting

Fast filling

### Graduated nebuliser cup guide


Faster preparation

### Reduced treatment times

May improve patient compliance



### HOT Top2 will operate in orientations from vertical to horizontal

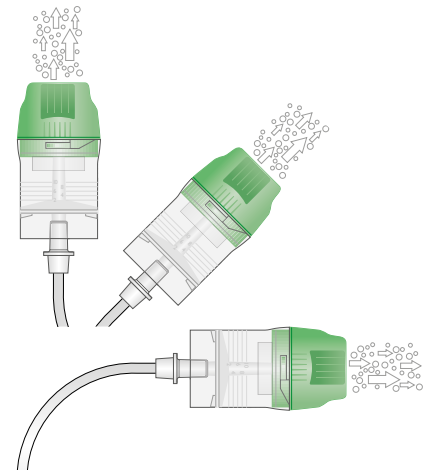

## Secure tube connection

### Sure-Loc™ base

HOT Top2 nebulisers are also available with a Sure-Loc base, when used with Intersurgical Sure-Loc tube & adaptor this system provides secure connections, reducing the risk of accidental tube disconnections from the nebuliser and the flow meter.

The oxygen tube screws onto the nebuliser base and the other end is connected to the flow meter with an adaptor.

The Sure-Loc tubing can also be used on a standard nebuliser spigot and Sure-Loc nebulisers will also accept standard tubing.



Push and turn to lock



## HOT Top™2 nebuliser, Intersurgical EcoLite™ mask kit

### Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

### Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patient's face.

### Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask body and seal, reducing the environmental impact.



### Intersurgical EcoLite

[Find out more](#)

### Eco range

[Find out more](#)



### Incurved nose seal

Conforms to different nose shapes designed to reduce the risk of aerosolised drug entering patient's eyes

### No metal nose clip

MRI compatible

### Soft face seals

Increased patient comfort

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

### A choice of "under chin" or "on chin" positions





Provides a better fit on a wider range of patient face shapes







## Product information and accesories

### HOT Top™2 nebulisers – Standard base

Code	Description	Tube length	Box Qty.
1472015	HOT Top2 nebuliser		65
1463015	HOT Top2 nebuliser, adult, Intersurgical EcoLite™ mask kit with tube	1.8 m	30 
1463016	HOT Top2 nebuliser, adult, Intersurgical EcoLite mask kit and tube with universal connector	1.8 m	30 
1468015	HOT Top2 nebuliser, paediatric, Intersurgical EcoLite mask kit with tube	1.8 m	36 
1468016	HOT Top2 nebuliser, paediatric, Intersurgical EcoLite mask kit and tube with universal connector	1.8 m	36 
1478015	HOT Top2 nebuliser mouthpiece kit with tube	1.8 m	45

### HOT Top2 nebulisers – Sure-Loc™ base

Code	Description	Tube length	Box Qty.
1471015	HOT Top2 nebuliser with Sure-Loc base		75
1473015	HOT Top2 nebuliser, adult, Intersurgical EcoLite mask kit with Sure-Loc base and tube	1.8 m	30 
1474015	HOT Top2 nebuliser, paediatric, Intersurgical EcoLite mask kit with Sure-Loc base and tube	1.8 m	36 
1475015	HOT Top2 nebuliser mouthpiece kit with Sure-Loc base and tube	1.8 m	45

Make an enquiry

 Lower environmental impact product

IS10.28 ▪ Issue 15 09.25



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

**Please think before you print**  
Save energy and paper.  
If you must print this information sheet please print it double sided.